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1. LEGAL INFORMATION

Operating instructions

Please read the operating instructions carefully before using the device for the first time. Toothbrushes may be used by children and persons with reduced physical, sensory or mental abilities or insufficient experience or knowledge, provided that they are supervised or instructed in the safe use of the appliance and understand the potential hazards. Cleaning and maintenance by the user must not be performed by unsupervised children. The device should be used only with the power supply or USB charger provided with the product and the brush sold by Y-Brush. Connect the charger only to an AC outlet (100-240V). The product should only be operated with the safety extra-low voltage corresponding to the marking on the product. Children should not play with the product.

IMPORTANT SAFEGUARDS

SAVE THESE INSTRUCTIONS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING

DANGER- To reduce the risk of electrocution:

1. Do not use while bathing.
2. Do not place or store product where it can fall or be pulled into a tub or sink.
3. Do not place in or drop into water or other liquid.
4. Do not reach for a product that has fallen into water. Unplug immediately.

WARNING

To reduce the risk of burns, electrocution, fire, or injury to persons :

1. Close supervision is necessary when this product is used by, on, or near children or invalids.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the product to a service center for examination and repair.
4. Keep the cord away from heated surfaces.
5. Never use while sleeping or drowsy.
6. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

INTRODUCTION

Y-Brush is a device designed and developed by oral health specialists to facilitate dental brushing for both children and adults. It allows you to brush your teeth in a simple, effective way, and in a few seconds thanks to a patented system.

This device is composed of a flexible brush that fits perfectly in the jaw. Thanks to the nylon fibers lining the brush and the vibrations of the handle, the teeth are brushed simultaneously and efficiently.

Remember that the recommendations in terms of oral health are to brush your teeth at least twice a day, in the morning and in the evening. An effective brushing helps to prevent the consequences related to a bad oral health. (periodontal, cardiovascular diseases, etc.) It helps to remove plaque, but is not intended to remove tartar deposits. It is recommended to see a dental surgeon every six months.

DIRECTION FOR USE

+ CHARGING DEVICE

- Before each use of the Y-Brush, make sure it is fully charged. In normal use, the indicator lights on the handle are lit in white. They will light in red when the battery level is too low and in green when the Y-Brush is plugged into the power supply.
- If the battery level of the handle is too low:
 - ♦ Connect the charger to a 110-220V power outlet or use the USB charger
 - ♦ Plug the charging cable into the part provided for this purpose on the handle. An illuminated halo around the on/off button lights up cyclically and indicates that the device is being charged. If the halo stays on green continuously, it indicates that the handle is fully charged. A period of 8 hours is sufficient to fully charge the handle. It will then have an autonomy of up to 3 months at a rate of 2 brushes per day.
 - ♦ The handle is ready to be used. It is imperative to close the cover that protects the connection of the recharging cable on the handle side to guarantee its watertightness.
 - ♦ The handle cannot be used during charging.

WARNING

THE BRUSH

- The appliance and its brush have the sole function of brushing teeth, they should not be used for anything else. Do not attempt to swallow all or part of the brush.
- The handle must not be inserted in the mouth without a brush.
- Make sure that the nylon fibers of the brush reach all areas of all teeth once in the mouth. If any area are not reached, brush them separately.

- Be aware of the risk of pinching when the brush is attached to the handle.
- The brushes are individual and have a maximum service life of 6 months under normal use. If there is any doubt as to whether the brush belongs to a third party and/or can be used by a third party, change the brush to avoid the risk of cross-contamination.
- Replace the brush after severe impact (even if there is no visible sign of damage).

+ POWER SUPPLY & CONNECTION

- Use the unit at room temperature (0-40°C). Avoid prolonged exposure of the product to sunlight or any other heat source. Do not use the charger near water (bathtub, washbasin, etc.) or outdoors. Never immerse or rinse the charger.
- Never touch the power cable with sharp or pointed objects. Never touch the plug with wet hands, and never unplug by pulling on the cable. Regularly check the cable for cracks or damage.
- Be sure to close the cover that protects the handle charging plug connection so that it remains watertight. The handle is only splash-proof.
- Do not place in a microwave oven or any other household appliance.
- After use, rinse the brush head thoroughly, shake off excess water and store it on the brush holder with the connector facing up.
- Do not dispose of the product in the household garbage. The handle has batteries and should therefore be taken to a collection point for electronic products. The charger must also be taken to such a collection point. Do not attempt to remove the batteries.

+ OTHER IMPORTANT INSTRUCTIONS

- **CAUTION:** Not suitable for children under 4 years old. The elements, in both the adult and child version, contains small parts that may cause a choking hazard. Therefore, they must not be used by a child without adult supervision.
- Make sure that the use of the device is not unpleasant and does not impair dental health. If in doubt, contact a dentist. If you have dental problems, wear braces, have dental implants or crowns, consult your dentist or physician before using the product.
- Do not open the product. Any warranty is void if the product is opened. Improper repair is dangerous for the user.
- Never use a product that has been opened or damaged.
- Do not wash the product in the dishwasher.
- People with pacemakers or other electronic implants should consult their physician beforehand.
- The product should not be used for commercial purposes.
- Keep the packaging out of the reach of children.
- The manufacturer is not liable for any failure to observe these operating instructions, for repairs that are not carried out in accordance with the state of the art, or for improper use of the device.

Attention,

Please note that the warranty does not apply if water is found in the handle or if the handle has not been used according to these instructions. It is therefore important to close the connector cover and not to leave the handle within reach of a sink, shower or bath.

- This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

- Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

- **Note:** This equipment has been tested and found to comply with the limits for a class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interferences to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna
- Increase the separation between the equipment and receiver
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.
- This digital apparatus does not exceed the Class B limits for radio noise emissions from digital apparatus set out in the interference-causing equipment standard entitled: "Digital Apparatus" ICES-003 of the Canadian Department of Communications.
- This device complies with Industry Canada License-exempt RSS standard(s).

Operation is subject to the following two conditions:

1. This device may not cause interference
2. This device must accept any interference, including interference that may cause undesired operation of the device.

FasTeesH SAS

14-22 Avenue Barthélemy Thimonnier, 69 300 Caluire-et-Cuire

E-mail: customers@y-brush.com

Tel: **+33783477158**

Site: www.Y-Brush.com

2. PACK CONTENT:

1 HANDLE / 1 BRUSH HEAD / 1 CABLE

and optional

1 STORAGE STAND / 1 STICKERS BOARD


1 TOOTHPASTE APPLICATOR (*child or adult*) /

1 TRAVEL BAG


3. SPECIFICATIONS

Discover your new toothbrush
with sonic vibrations.

BRUSH HEAD

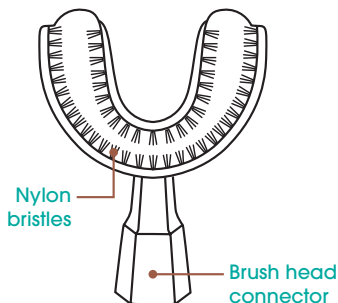HANDLE

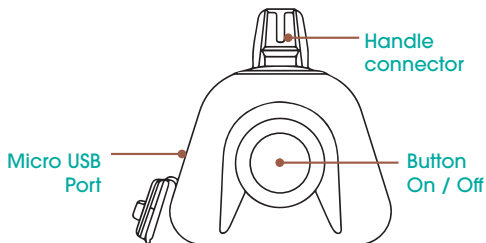Splash-proof handle.
See p12

YOU WILL NEED AN ADAPTATION TIME OF 1 TO 2 WEEKS. FOLLOW THE GUIDE!

4. GET STARTED WITH YOUR NYLONSTART TOOTHBRUSH

The Y-Brush NylonStart toothbrush allows you to adjust the brushing time according to your needs.

MODE

• Press the center button to start or stop your brushing.

WEEK 1

For this first week, **take the time** to chew well and to make the necessary number of oscillations to have an effective brushing. Your brushing time will not exceed 30 seconds. Please refer to **chapter 6** for the Instructions for use.

WEEK 2

Once you have understood and adopted the technique, you can **speed up the movement**.

WEEK 3

Now that you are an expert, you can adopt a more tonic gesture.

Don't hesitate to gauge the speed of the oscillations to ensure smooth teeth after each brushing.

⊕ = MODE INTENSITY

5. GET STARTED WITH YOUR NYLONBLACK TOOTHBRUSH

This range offers you **4 MODES**: **soft** mode, **normal** mode, **intensive** mode and **unlimited** mode.

4 MODES

- **SHORT PRESS** the center button to select your mode.
- **LONG PRESS** the center button to start or stop your brushing.

WEEK 1

Start with the **unlimited mode**. It allows you to take the time to learn the right movements, while ensuring effective brushing.

The process remains the same : **5 BITES** and **5 OSCILLATIONS** for one jaw.
Be patient and you'll see that the results will follow.

WEEK 2

Switch to **soft** mode, which lasts **15 seconds**.
Do not hesitate to speed up the movement.

WEEK 3

You can now switch to the **normal** mode, which is the 10-second mode. The oscillations should be faster. Your brushing will be more tonic and even more efficient.

IN A HURRY? The **NYLONBLACK** toothbrush has a **5-second intensive mode** to be used with the same number of oscillations and bites.

6. INSTRUCTIONS FOR USE


3 steps to a fast and effective brushing

6.1 APPLY YOUR TOOTHPASTE

With the **Y-BRUSH APPLICATOR**: place it on the toothpaste. Spread a dab of toothpaste along the entire length of the brush.

On the contrary, if you don't have any, put a dab of toothpaste on your tongue and spread it over all your teeth.

After applying the toothpaste, remember to **WET** the brush head.

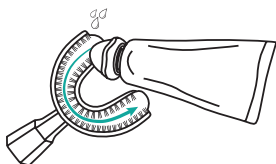TOOTHPASTE TABLET

Crunch the tablet and spread the pieces over your teeth.

6.2 THREE, TWO, ONE... BRUSH!

In our user guide, we talk about **OSCILLATIONS**. An oscillation is a **BACK AND FORTH MOVEMENT** with your toothbrush. This movement is to be done while brushing, oscillating from left to right, then from right to left, repeatedly.

2 RULES TO REMEMBER FOR A BRUSHING CYCLE

5

&

10

OSCILLATIONS


BITES

PER OSCILLATION

- 1 **SELECT YOUR MODE** (ONLY FOR THE **NYLONBLACK** VERSION)
- 2 **PLACE THE BRUSH ON ONE OF YOUR JAWS**
- 3 **PRESS THE START BUTTON TO START THE VIBRATIONS***
- 4 **BITE AND OSCILLATE THE BRUSH**
- 5 **DO THE SAME FOR THE OPPOSITE JAW**

RUN YOUR TONGUE OVER YOUR TEETH: IF THEY ARE SMOOTH, THE PLAQUE HAS BEEN REMOVED.

**REFER TO CHAPTER 4 OR 5 CORRESPONDING TO YOUR PRODUCT*

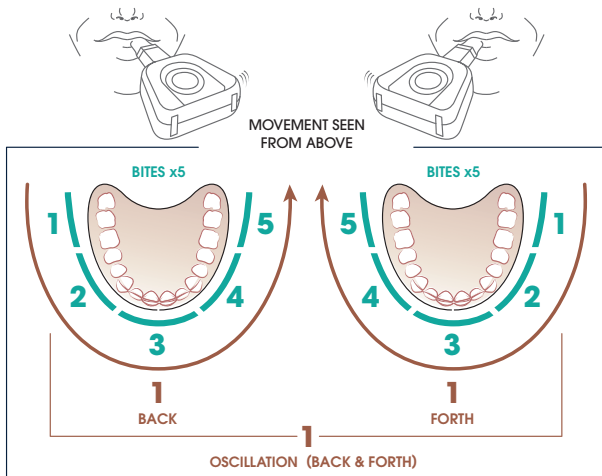
Explanatory scheme of **AN OSCILLATION**, that's to say a back and forth movement from left to right and then from right to left, accompanied by 10 bites.

6.3 CLEANING AND STORAGE

When you have finished brushing your teeth, pinch the base of the brush to unlock it, then **RINSE** it with clean water. Run a finger through the brush to remove the toothpaste. Let it dry on its **STAND** (optional), or upside down on the edge of a glass.

To take care of your brush head, we recommend that you soak it once a week in a diluted mouthwash.

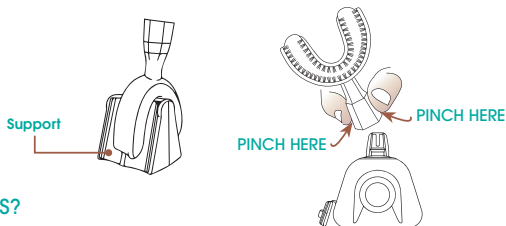ANY QUESTIONS?

Write us at:

customers@y-brush.com

Our customer service is committed to respond within 24 hours to give you the best experience possible.

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DISCOVER

our video tutorials by scanning this QR code,
or visit y-brush.com/support

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